

Barella Shower Trolley Range

The feature
packed mobile
shower trolley



Welcome to Wealden Rehab



Simon Campbell
Managing Director

Over the years we have dealt with many professionals, working with users of showering and bathing equipment across the UK. One thing we have heard many times over is the importance of shower trolleys that are robust and carer friendly – they should be comfortable and have a good platform size to enable turning and aid personal hygiene routines.

The Barella trolley range has been designed to combine all of these important features, and provides an ever-expanding collection of trolley sizes and options to bring you the perfect solution for shower trolley bathing. Many customers have taken advantage of our quick delivery service on the most popular models in the range, from high stock levels in our central warehouse.

The Barella range is also available through the NHS Supply Chain - we know and can demonstrate to our customers that our products will perform in the most demanding of environments.

Whatever your requirements, contact one of our product advisors. You can be sure that when you partner with Wealden Rehab you will get a solution to be proud of. We look forward to hearing from you!



Care Equipment Excellence

The Barella Shower Trolley Range

Welcome to the world of Barella,
where your comfort is our focus,
and your investment is our passion.

Barella is the product of many years of research, fine tuning and critical observation. The latest range of shower trolleys is designed to offer stability, comfort and convenience.

Developed into a modular selection of chassis, platforms and options, your choice of product is both extensive and focussed. We now believe we have the broadest collection of shower trolleys to choose from, ensuring that 'making do' is a sentiment of the past.

At the heart of the Barella range is the ability to choose the specification of your trolley according to your individual requirements.

Modular accessories can be chosen at the time of order, making the breadth of choice the widest in the market.




Barella Platform Size

Paediatric

Inside width 60 cm


For use especially with children up to 160 cm tall.

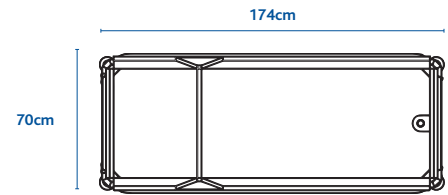

Adult

Inside width 68 cm

For users up to 188cm tall.

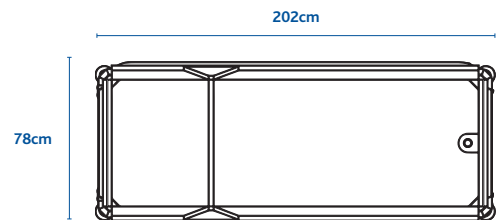


Bariatric

Inside width 80 cm

For larger users up to 188 cm tall.

*This wide platform is also available on the Adult Barella

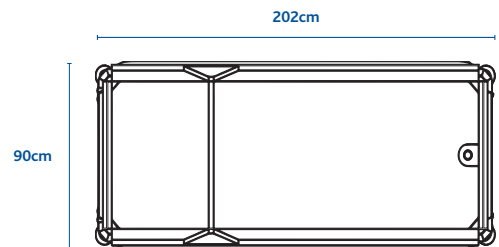

Height Adjustment

Adjusting the shower trolley platform

The shower trolley height can be adjusted via electric actuator or hydraulic piston.

The electric version can be adjusted using a pedal control or a handset.

Electric versions are supplied with a battery charger




Hydraulic, with double pedal

Adjustment via hydraulic piston, with pedals on both sides of the trolley.



Electric, with handset

In this case, the carer can adjust the height electrically using a handset.

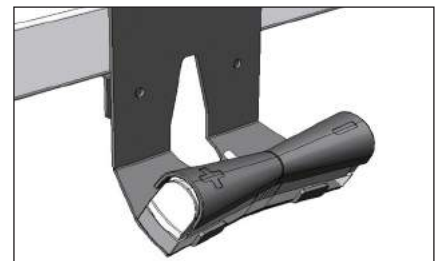

Electric, with pedal

This is an alternative option for electric height adjustment. The control is securely placed in a support, which can be attached at the base of the structure if so required by the user. The control can also be used as handset, or removed from its support and attached wherever needed using the magnet.

Angle Adjustable Backrest

Angle-adjustable backrest

This helps give greater comfort to both the user and carer during washing, as it allows to change the patient's position. The backrest can be easily raised by the carer via gas spring and ergonomic control.

A. Non-adjustable backrest

The patient remains in a supine (flat) position inside the mattress. The patient's head can be lifted using a wedge pillow, available as an accessory.

B. Adjustable backrest

The backrest can be raised to a 63° angle. The backrest can be locked in any intermediate position. A separate accessory can also be requested to provide greater head support.



Extended Head Support

An optional extension will support the users head when the backrest winto position once the backrest has cleared the endrail.

Castor Choice

A. Castors with 12.5cm diameter, individually braked

Castors have individual brakes, and one castor also has a "directional lock" function.



Castor with individual brake



Castor with directional lock

B. Castors with 15cm diameter, central braking system

These castors have a larger diameter with a central braking system and directional lock function. These functions can be operated from either side of the trolley. one castor also has a "directional lock" function.



Central braking activated



Directional lock activated

Accessories



BH733 – Spare battery

This is a rechargeable, 24V spare battery recommended to prevent any interruptions during use.



B2WP – Wedge pillow

Pillow with 15° tilt and dimensions cm 55 x 55 x 14½ h.



B2SCP – Semi-cylindrical pillow

Dimensions: 55 x 15 x 7 cm

Barella Paediatric

A feature packed
mobile shower trolley



Based around a smaller chassis, and a reduced weight capacity, the Paediatric model is aimed to accommodate children and small adults within a compact footprint.

The Paediatric trolley is available with either hydraulic or electric adjustment, and can be equipped with an elevating backrest.

Technical Information

Height adjustment Control	Hydraulic/Battery
Maximum Safe Working Load	170kg
Trolley Weight	80kg
Castor diameter	12.5cm
Castor Locking System	Individual
Directional Locking Castor	Standard
Folding side rails	Standard
Trendelenburg Platform Tilt to 12°	Standard
Water Drainage Platform Tilt to 5°	Standard
Drain Hose Length	150cm
IP Protection	X5
Power supply	220v, 500Ma to battery charger
Electrical Safety	Onboard circuit - 24v

Order Codes

With rising backrest

Electric

BP31C1121

Hydraulic

BP31D1121


Accessories

- BH733 Spare battery pack, 24v, rechargeable
- BH744 Battery charger, wall-mounted
- B2WP Wedge Pillow
- B2SCP Semi-cylindrical Pillow
- B2RH Replacement Handset



The backrest tilt can be adjusted up to a 63° angle via ergonomic control.

Key Dimensions (in mm)

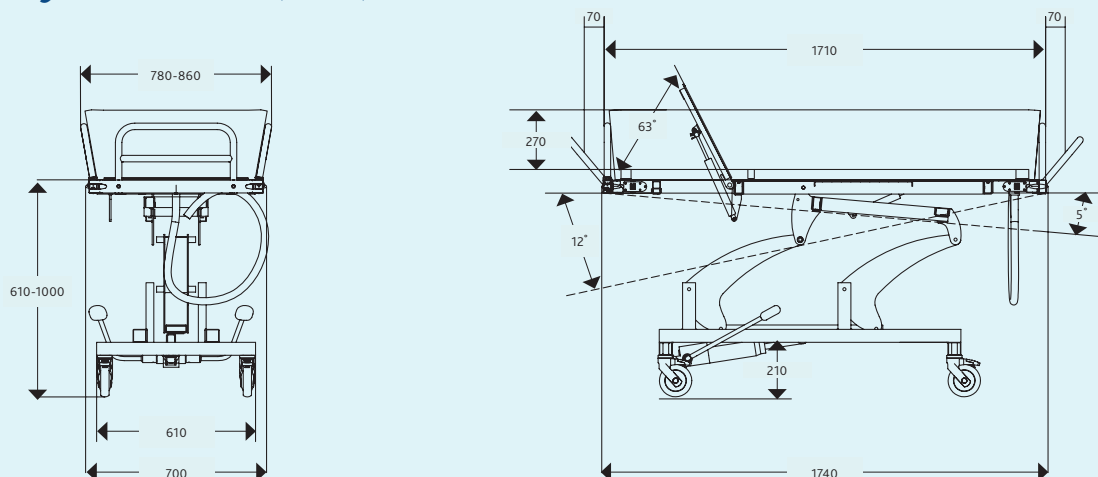

Specification Choices

HEIGHT ADJUSTMENT

- Hydraulic, with double pedal
- Electric, with handset
- Electric, with pedal

PLATFORM

- Inside width 60cm

LENGTH EXTENSION

- Fixed rails on head and foot end
- Head end rail folding
- Rails on head and foot end folding

BACKREST ADJUSTMENT

- Flat Platform
- Rising Backrest

BRAKING SYSTEM

- Individual on all castors, 12.5cm diameter

Barella Adult

A feature packed
mobile shower trolley



Our most popular model, the Adult chassis has a superb lifting frame which guarantees stability at every working height. Many options are available on this model, making it the trolley of choice for demanding environments.

An increased weight capacity of 180kg makes this an excellent trolley for everyday showering. Choose central locking castors for increased ease of use.

Technical Information

Height adjustment Control	Hydraulic/Battery
Maximum Safe Working Load	200kg
Trolley Weight	106
Castor diameter	12.5cm / 15cm
Castor Locking System	Individual / Central
Directional Locking Castor	Standard
Folding side rails	Standard
Trendelenburg Platform Tilt to 12°	Standard
Water Drainage Platform Tilt to 5°	Standard
Drain Hose Length	150cm
IP Protection	X5
Power supply	220v, 500Ma to battery charger
Electrical Safety	Onboard circuit - 24v

Order Codes

Standard	BR32D2111	BR32C2111
Cent locking/rising backrest	BR32D2122	BR32C2122
Wide	BR32D3122	BR32C3122



Accessories

- | | |
|-------|---------------------------------------|
| BH733 | Spare battery pack, 24v, rechargeable |
| BH744 | Battery charger, wall-mounted |
| B2WP | Wedge Pillow |
| B2SCP | Semi-cylindrical Pillow |
| B2RH | Replacement Handset |



The rail width can be reduced to 78 cm.

Key Dimensions (in mm)

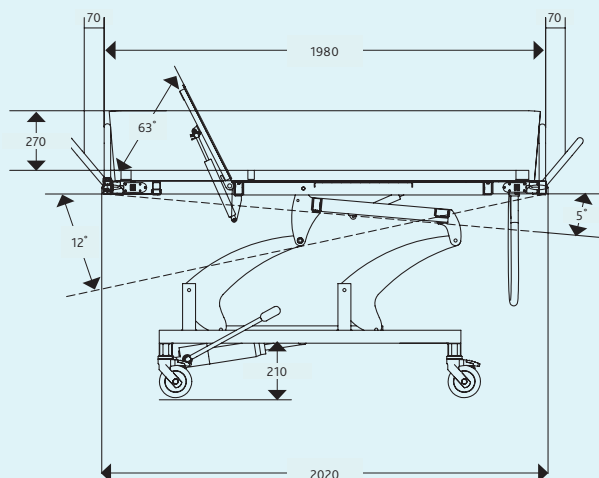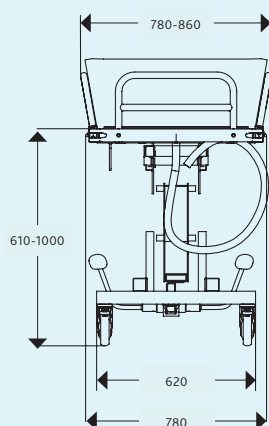

Specification Choices

HEIGHT ADJUSTMENT

- Hydraulic, with double pedal
- Electric, with handset
- Electric, with pedal

PLATFORM

- Inside width 68 cm
- Inside width 80cm

LENGTH EXTENSION

- Fixed rails on head and foot end
- Head end rail folding
- Rails on head and foot end folding

BACKREST ADJUSTMENT

- Flat Platform
- Rising Backrest

BRAKING SYSTEM

- Individual on all castors, 12.5cm diameter
- Central-braking with castors, 15cm diameter

Barella Bariatric

A heavy duty
feature packed
shower trolley



For the more substantial user, the Heavy Duty chassis delivers real performance. Twin electric actuators combine to provide effortless lifting and platform tilt. A detachable battery pack provides power whilst the trolley is moved around.

A wide platform top gives adequate space for hygiene and handling, and the adjustable platform tilt provides extra comfort for the user enabling safe transfers, but also greater functional independence.

Technical Information

Height adjustment Control	Battery
Maximum Safe Working Load	225kg
Trolley Weight	106
Castor diameter	15cm
Castor Locking System	Central
Directional Locking Castor	Standard
Folding side rails	Standard
Trendelenburg Platform Tilt to 12°	Standard
Water Drainage Platform Tilt to 5°	Standard
Drain Hose Length	150cm
IP Protection	X5
Power supply	220v, 500Ma to battery charger
Electrical Safety	Onboard circuit - 24v

Order Code

BB33



Accessories

- BH733 Spare battery pack, 24v, rechargeable
- BH744 Battery charger, wall-mounted
- B2WP Wedge Pillow
- B2SCP Semi-cylindrical Pillow
- B2RH Replacement Handset



Rails can be dropped with one hand.
Mattress height is 27 cm.

Key Dimensions (in mm)

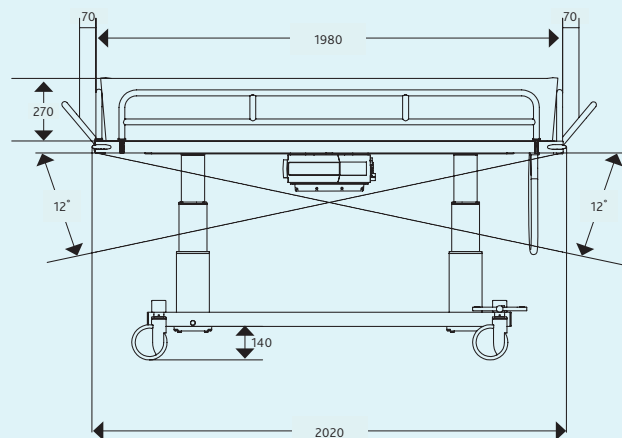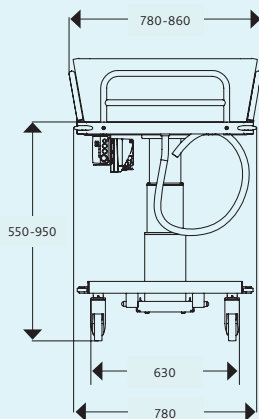

Specification Choices

HEIGHT ADJUSTMENT

- Electric, with handset

PLATFORM

- Inside width 68 cm
- Inside width 80cm

LENGTH EXTENSION

- Fixed rails on head and foot end
- Head end rail folding
- Rails on head and foot end folding

BACKREST ADJUSTMENT

- Non-tilting

BRAKING SYSTEM

- Central-braking with castors, 15cm diameter



0845 658 8411
hello@wealdenrehab.com
www.wealdenrehab.com



113 Hopewell Drive, Chatham, Kent ME5 7NP | Wealden Rehab is a trading name of Wealden Rehab Ltd

