



# Heron

A stainless steel pool hoist which can be configured for different environments



# Welcome to Wealden Rehab

Thank you for your interest in Wealden Rehab. Since 2003 we have become established as a leader in the supply of specialist care equipment solutions.

Wealden Rehab is a family managed business with a passion for service that will make a real difference to the lives of our users and carers everywhere. We work hard to ensure our products are delivered with the best possible experience.

## Our mission statement

'Superior solutions and innovative design which will improve the quality of care for our clients, underpinned with simple honesty, enthusiastic service and integrity at every step.....'

Our portfolio covers products in bathing and hygiene, patient handling, and postural and pressure management. We pledge to offer specific assessments, rapid itemised quotations, and a personal delivery service.

We have long experience in working with the NHS, healthcare professionals (occupational therapists in NHS, social services and private practice), architects, builders, and nursing and residential homes as well as the friends and families of those that need our products. Whoever you are we look forward to assisting you in finding the right solution.

# Heron

Heron is designed to provide a safe and secure transfer into a pool. Heron is available in four different versions, allowing for pools built above or below the ground, and the user can be transferred in a fabric sling or on a solid seat.

The Heron is constructed from stainless steel which enables it to be confidently installed in harsh environments. The mast of the hoist can be rotated through 360° enabling a user to be easily hoisted from poolside, and placed in a convenient position within the pool.

The Heron uses a 24V detachable and rechargeable battery system to provide smooth and trouble-free lifting power for users up to 135kg/21st.

A choice of two fixing methods ensures that the Heron can be deployed in the majority of situations. A surface mounting plate is ideal where there are underground services, whilst a drilled socket is perfect for maintaining clean poolside design.



## DESIGNED FOR POOLS

Above or below ground, with or without barriers.

## STAINLESS STEEL CONSTRUCTION

For reliability and durability in all environments

## CHOICE OF PATIENT INTERFACE

Choose a sling spreader bar or detachable patient seat

## BATTERY POWERED LIFT

For effortless operation

## CHOICE OF MOUNTING PLATES

Depending on the pool design and surroundings



### POOL SOCKET


Concealed fixing consisting of a stainless steel socket, which provides a secure housing for the pivot pin. Complete with cap.



### POOL PLATE

Surface mounted fixing method by means of a stainless steel plate. An ABS cover is provided.

## DIMENSIONS

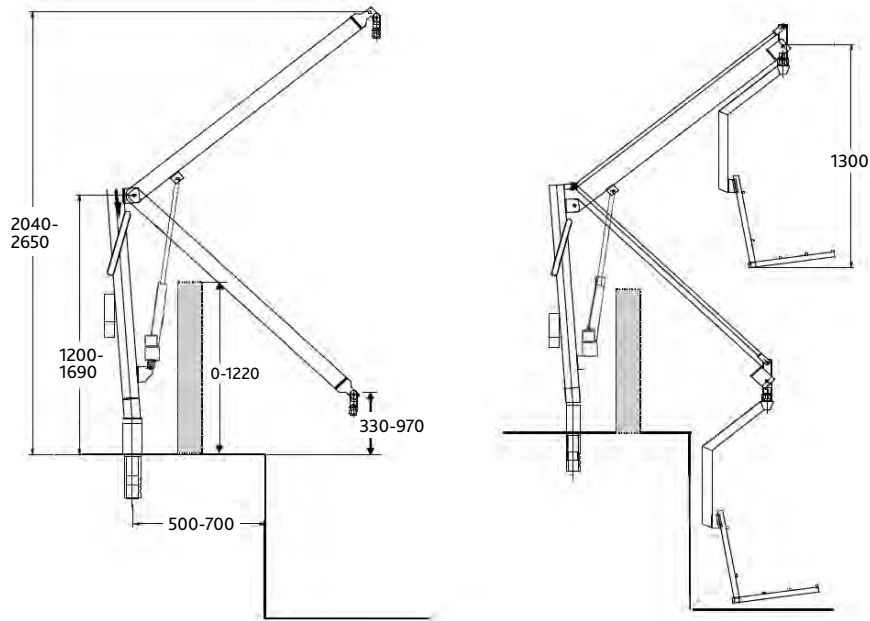

Max user weight	134kg / 21st
Battery capacity	24v

Dimensions in mm

For more information or for an assessment  
call 01634 813388 or email [sales@wealdenrehab.com](mailto:sales@wealdenrehab.com)  
Visit [www.wealdenrehab.com/heron](http://www.wealdenrehab.com/heron)



## Innovative care equipment with professional installation

113 Hopewell Drive, Chatham, Kent ME5 7NP | Wealden Rehab is a trading name of Wealden Rehab Ltd

Our brands

**Barella™**  
mobile shower trolley

**Valencia**  
specialist riser chair

**boss™**  
uniquely shaped for support

**RAZ**  
modular shower chair

**sedoccia**  
specialist hygiene chair

**KINETIC BALANCE**  
postural positioning belts

**AltoSlings™**  
support you can trust

**heron™**  
professional pool hoist